


HOSPITAL BEDS



27841


Height adjustable
160 to 245 cm

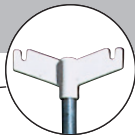

27610 + accessories

Lifting pole
Code 27644
(not included)

Collapsible support
code 27641
(not included)

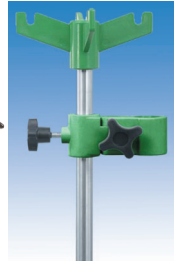

27646

BED WITH ADJUSTABLE HEAD SECTION

- 27600 BED WITH 1 ARTICULATION
- 27601 BED WITH 1 ARTICULATION - with wheels

Painted steel frame with net.

1 Pipe frame with net for head movement, 200x85 cm.

2 Heads in chromed pipe complete with MDF bilaminated coated panel. Movement by a crank. Protection of end run.

Rubber wheels of Ø 100 mm, with brakes (only for code 27601)

Size: 206x89xh 90 cm. Height of frame from bottom 52/55 cm.

Size of frame: 200x85 cm. Weight 32 kg.

Made in Italy

BED WITH ADJUSTABLE HEAD AND FOOT SECTION

- 27610 BED WITH 3 ARTICULATIONS
- 27611 BED WITH 3 ARTICULATIONS - with wheels

Painted steel frame with net.

3 bed frames with nets in painted steel pipe, 200x85 cm.

2 Heads in chromed pipe complete with MDF bilaminated coated panel. Movement by 2 cranks. Protection of end run.

Rubber wheels of Ø 100 mm, two of which with brake (only 27611)

Size 206x89xh 90 cm. Height of frame from bottom 52/55 cm.

Size of frame: 200x85 cm.

Weight 42 kg. Made in Italy

- 27605 STANDARD BED - no castors
- 27606 STANDARD BED - Ø 50 mm castors

Simple bed made of epoxy - coated steel. Head and foot ends with plastic feet. Rigid mesh mattress platform.

- 27607 STANDARD PLUS BED - no castors
- 27608 STANDARD PLUS BED - Ø 100 mm castors

Chrome-plated steel bed ends. Laminated plastic panels.

Rigid mesh mattress platform. Adjustable back-rest

Size: 206x89xh 90 cm. Made in Italy.

- 27841 I.V. STAND - 5 wheels trolley



27605



27608

GIMA code	ACCESSORIES FOR ALL BEDS
27641	Collapsible supports (for all types of bed)
27644	Lifting pole
27646	4 Hook IV stand for bed - stainless steel

BEDSIDE TABLES - MADE IN ITALY



27665

Size: 42x40xh 83 cm



27669

Size: 42x40xh 83 cm



27674

Size: 42x40xh 80 cm



27662

Size: 42x52.5xh 80 cm



27679

Four revolving castors Ø 60 mm (2 with brake) with bumpers Ø 100 mm.



27678

DOUBLE FACE

BEDSIDE TABLES

Epoxy powder coated steel sheet tube or body.

Thermoformed ABS top and drawer.

Laminated plastic drawer face and door (27665 and 27666 only) 10 mm thick.

Bedside table with door has inner removable shelf and is available in left (L) or right (R) version.

- 27669 BEDSIDE TABLE WITH DRAWER AND LEFT DOOR
- 27665 BEDSIDE TABLE WITH DRAWER AND LEFT DOOR
- 27666 BEDSIDE TABLE WITH DRAWER AND RIGHT DOOR

FEBO BEDSIDE TABLES

Epoxy powder coated steel sheet 8/10 body. Towel holder and lateral bottle holder (27662 only). Soft rubber antitrauma handle. Sliding drawer with removable ABS tray. Inner removable shelf. 4 revolving castors (2 with brake)

- 27674 FEBO BEDSIDE TABLE - blue avio RAL 5014
- 27675 FEBO BEDSIDE TABLE - other colours*

RAL 5010	RAL 2003	RAL 6019	RAL 1013	RAL 7040
----------	----------	----------	----------	----------

*Available in 30 days. Minimum order 10 pcs.

- 27662 FEBO BEDSIDE TABLE WITH BOTTLE HOLDER - blue avio RAL 5014

BEDSIDE TABLE WITH DRAWER, DOOR, BOTTLE - GLASS HOLDER, OVERBED TABLE

Epoxy powder coated steel sheet body. Thermoformed ABS top with towel holder, lateral bottle and glass holder. Sliding drawer with removable ABS tray. Laminated plastic drawer face and door 10 mm thick. Anti-trauma handles.

Revolving, extensible, adjustable height cm 90/110 ABS overbed table top. Removable inner shelf. Plastic grid at bottom.

- 27678 BEDSIDE TABLE - DOUBLE FACE
- 27679 BEDSIDE TABLE - SINGLE FACE