

DONOR ARMCHAIR

- 27557 DONOR ARMCHAIR - blue - mechanical
- 27563 DONOR ARMCHAIR - blue - electrical
- 27567 DONOR ARMCHAIR WITH WHEELS - white

Strong and stylish armchair with structure made with oval, round and square metal section tubing (legs oval tubing 50x25x1.5; central frame square tubing 30x30x1.5, cushions round tubing 20x1.5 diam.). Structure fully coated with white epoxy powder. High density polyurethane foam platform with wooden base covered by leatherette upholstery.

Both materials are flame resistant Class 1 IM. Synchronized adjustment of the two cushions by a gas spring, cylinder/piston, (27557, 27567) or by a linear actuator (27563) with a push-button panel, allow a quick moving into TR/RTR position. The backrest inclination ranges from +60° to -5° in Trendelenburg position, leg rest from -80° to +5°.

The arm-rests can be height adjusted and are reclinable. When removed, they can be rotated up to 180°.

Four adjustable feet ensure maximum stability (27557 and 27563) or 4 castors, 2 with directional brake (27567) ensure mobility. Load: 150 kg. Size: 175x55xh 51 cm. Made in Italy.

- 27553 IV STAND
- 27556 COUCH ROLL HOLDER



27557



27557

White ral colour 9003 structure, minimum thickness 70/80 micron. Glossy finish.



WITH WHEELS



27567



27563

ELECTRICAL

DONOR WHEELCHAIR

- 27561 DONOR WHEELCHAIR

Donor wheelchair with structure made of epoxy painted oval and round steel tube.

Backrest adjustable by gas pump.

Padded leg rest manually adjustable.

Metal footrest extractable. Armrest with articulations, manually adjustable in height.

Padding made of high-density polyurethane foam, fireproof Class 1 IM. Covering made of washable, plastic material (leatherette), fireproof Class 1 IM.

Back push handle.

Mattress overall size: 166x55 cm. Carton size: 185x67xh 58 cm. Four wheels Ø 100 mm, 2 directional and 2 with brake.

Net weight: 53 kg. Gross weight: 55 kg. Maximum load: 150 kg. Supplied with blanket tray with belts to fix oxygen bottle, IV stand and chart holder support. Made in Italy.

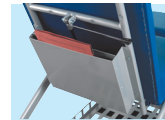


Chart holder support



27561

Blanket tray with belts to fix oxygen bottle (bottle not included)

HEIGHT ADJUSTABLE GYNAECOLOGICAL BED

- 27506 HEIGHT ADJUSTABLE GYN BED - green
- 27507 HEIGHT ADJUSTABLE GYN BED - blue
- 27496 HEIGHT ADJUSTABLE GYN BED - any colour*

*Available in 20 days, see below available colours

Strong and solid structure in oval tube 50x25x2 mm painted with epoxy powder, white colour.

The three sections top has high density padding, covered with washable, fire proof, leatherette material class 1 IM.


Full adjustable

Height, up to 95 cm, is adjustable by electric engine (6000 N power) with foot control, head section is adjusted by a gas cylinder (400 N power), central section by a gas cylinder (600 N power) up to 65°. This gyn bed is supplied with 2 quality foam leg holders, stainless steel tray and couch roll holder.

Maximum load: 150 kg. Weight: 70 kg. Made in Italy.

- 27509 QUALITY LEG HOLDER - couple - spare

Kit leg holder includes 2 leg holder and arms. Suitable for 27496, 27505, 27506, 27507.

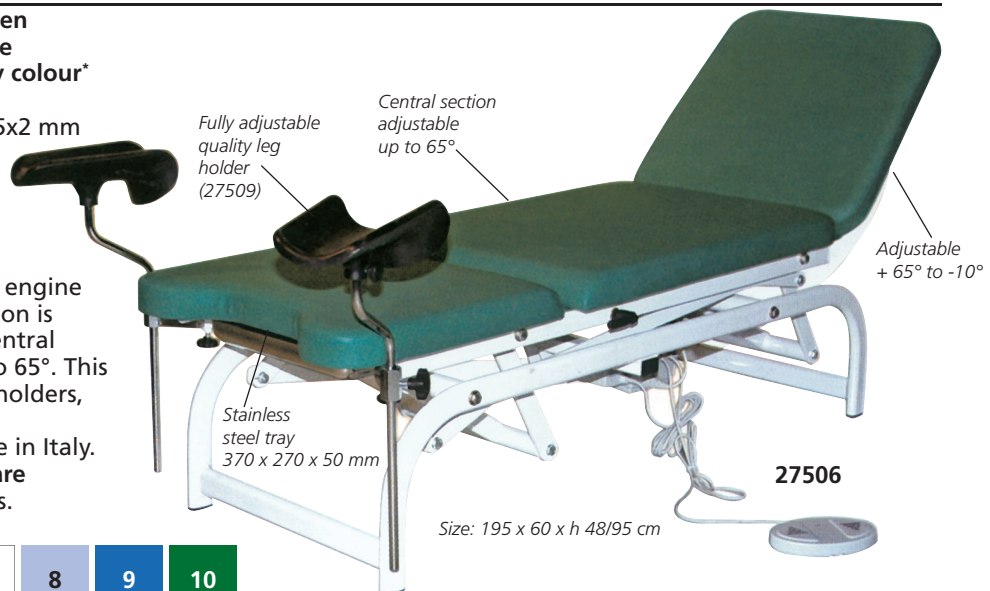

27506

Fully adjustable quality leg holder (27509)

Central section adjustable up to 65°

Adjustable + 65° to -10°

Stainless steel tray 370 x 270 x 50 mm

Size: 195 x 60 x h 48/95 cm

